2pro®


Total Access Disposable Propphy Angle

Disposable Propphy Angle

- FR Angle à prophylaxie dentaire jetable
- ES Ángulo profiláctico desechable
- DE Einweg-Prophylaxe-Winkelstück
- IT Innesto a contrangolo per profilassi monouso
- KR 일회용 프로피 앵글 (Disposable Propphy Angle)

 Manufacturer: Premier® Dental Products Company
1710 Romano Drive, Plymouth Meeting, PA 19462 U.S.A.
610-239-6000 • 888-670-6100
E-mail: info@premierdentalco.com • premierdentalco.com

U.S. Patent #7,104,796
U.S. Patent #7,153,133

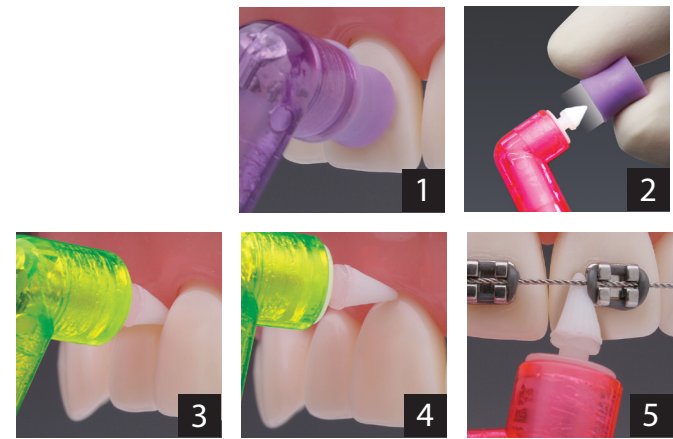EC REP MDSS GmbH
Schiffgraben 41
30179 Hannover, Germany

DFU# 55001 1219027 Rev10 JC

US / UK

2pro® Total Access Disposable Propphy Angle

Description:

2pro® is a single-use disposable propphy angle uniquely designed with a propphy cup and Total Access Tip for cleaning and polishing teeth during a dental hygiene appointment. 2pro's propphy cup is ideal for cleaning and polishing facial and lingual tooth surfaces. Remove the propphy cup and use the Total Access Tip to clean and polish gingival and interproximal areas, veneer margins and around implants or orthodontic brackets. 2pro's Total Access Tip can be used to clean and polish occlusal surfaces. 2pro can be used with Enamel Pro® Propphy Paste or any other propphy paste. 2pro is latex free and fits a nose-cone low-speed handpiece.

Indications for Use:

2pro is a single-use disposable propphy angle uniquely designed with a propphy cup and Total Access Tip for cleaning and polishing teeth. 2pro is designed for dental use only.

Contraindications:

There are no known contraindications for 2pro Total Access Disposable Propphy Angles.

WARNINGS:

2pro is designed for use on a low-speed handpiece. Do not operate 2pro above 3,000 rpm. **Do not reattach the cup after removing it from the hub of your 2pro Total Access Disposable Propphy Angle.**

PRECAUTIONS:

2pro is designed as a single-use, disposable device. Do not sterilize or reuse. Dispose of 2pro properly after each patient use.

Operating Instructions:

1. Remove 2pro Total Access Disposable Propphy Angle from packaging.
2. Insert the shaft of 2pro into the nose-cone of your low-speed handpiece. Align the slot of the propphy angle with the guide pin on your handpiece and lock into place.
3. Safely operate 2pro for a few seconds outside the mouth to ensure that it is attached securely to the handpiece.
4. To polish facial and lingual tooth surfaces, use the cup on your 2pro Total Access Disposable Propphy Angle. **(See Image 1)**
5. To clean and polish gingival and interproximal areas, veneer margins and around implants or orthodontic brackets, remove the propphy cup and use the Total Access Tip. **(See Image 2, 3, 4 & 5)**
6. To clean and polish occlusal surfaces, use the cup or Total Access Tip on your 2pro Total Access Disposable Propphy Angle.
7. Do not operate 2pro above 3,000 rpm.
8. **Do not reattach the cup after removing it from the hub of your 2pro Total Access Disposable Propphy Angle.**
9. 2pro is designed for single-patient use. Do not sterilize or reuse.
10. Dispose of 2pro properly after each patient use.

REF 2pro® with Total Access Tip - Pkg 144
5500101 – 2pro® Soft/Short Cup (Purple Cup)
5500102 – 2pro® Firm/Short Cup (Blue Cup)
5500103 – 2pro® Soft/Long Cup (White Cup)
5500104 – 2pro® Firm/Long Cup (Green Cup)

REF 2pro® with Total Access Tip - Pkg 1,000
5510001 – 2pro® Soft/Short Cup (Purple Cup)
5510002 – 2pro® Firm/Short Cup (Blue Cup)
5510003 – 2pro® Soft/Long Cup (White Cup)
5510004 – 2pro® Firm/Long Cup (Green Cup)

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