Premier® Universal Primer

Zirconia and Glass Ceramics, Metal and Composite Primer

Directions For Use

FR Mode d'emploi
ES Indicaciones de uso
DE Gebrauchsanleitung
IT Istruzioni per l'uso

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
Manufacturer:
Premier® Dental Products Company
1710 Romano Drive, Plymouth Meeting, PA 19462 U.S.A.
(888) 670-6100 • (610) 239-6000 • premierdentalco.com

To obtain an SDS visit premierdentalco.com or call Premier at 610-239-6000.

Made in Israel.

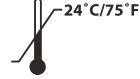Rx ONLY

DFU#31561 • 1119058 Rev1 BJM

US / UK

Premier® Universal Primer

Zirconia and Glass Ceramics, Metal and Composite Primer

Product Description:

Premier Universal Primer is a single component universal primer that provides an enhanced adhesive surface to zirconia, glass ceramics, metal and composite materials. The primer does not require additional conditioning such as sandblasting or acid etching. Its low viscosity and excellent wetting properties promote ideal surface coverage.

Indications:

Premier Universal Primer is a prescription-only, one-step silane coupling agent intended for preparation of ceramic material such as porcelain, ceramic, hybrid ceramic or composite and zirconia restorations designed to improve adhesion to resin material prior to cementation or restoration. This primer does not require prior sandblasting or the use of an etchant.

Precautions:

- Uncured methacrylate resin may cause contact dermatitis and damage the pulp. Avoid contact with skin, eyes and soft tissue. Wash thoroughly with water after contact.

- For intraoral repairs properly isolate the restorative area with a rubber dam.
- Avoid contact with the skin; unpolymerized (meth)acrylate resins may cause skin sensitization in susceptible persons. In case of contact, wash skin with soap and water.
- Keep out of the reach of children

Contraindications:

Do not use on patients with known sensitivity or allergy to methacrylate resins. Wash exposed area completely to remove material. If an allergic reaction is observed, consult a physician immediately and discontinue use with this patient.

Note:

Sandblasting the bonding surface of zirconia, metal and composite restorations prior to using Premier Universal Primer can be used to help achieve maximum bond strength for cases that require additional retention (e.g., short and/or over-tapered preparations).

Instructions For Use:

- Clean the internal surface of the restorations with alcohol and dry well with oil free air.
- Apply a generous amount of Premier Universal Primer onto the

- surface to be repaired with an applicator brush.
- Wait one minute.
- Dry gently for 20 seconds with oil-free air.
- Cement the restoration using a resin-based definitive cement according to manufacturer's instructions for use.

Storage:

Keep away from direct sunlight. No refrigeration required. If refrigerated, allow the bottle to reach room temperature before use. Do not freeze. Do not use after expiration date. See expiry date on the label. Shelf life is 2 years after the date of manufacturing.

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