

WOW!

Self-articulated and precise models in two fast pours!



W.O.W.[™] & W.O.W.²

The W.O.W. and W.O.W.² are innovative single-use articulators. Now you can articulate the case exactly as the dentist took the impression. Use for all types of cases, dual impression trays, pre-poured models and conventional impressions. The W.O.W. or W.O.W.² can provide a self-articulated, dependable and sturdy model in two fast pours!

Save money and time!

- No need to track and clean metal articulators
- No accelerator or glue required
- Two pours per case
- Minimizes grinding and trimming

Intelligent design provides accurate dies and models

- Clear plastic view allows for accurate pin placement
- Pin placement design eliminates rotation
- Shallow retention grooves minimize dimensional change from stone expansion
- Fixed hinge allows for positive occlusal relationship between upper and lower models and appropriate lateral movement
- Remount bite, if necessary, with heat in less than a minute

W.O.W.² features Flash Air Dam[™] technology

- Keeps unwanted debris from lodging in pin holes
- Keeps cases neat and delivers professional presentation



Available in full arch and posterior



Fast, easy technique for
accurate dies and models



W.O.W.™ & W.O.W.2




1. Trim excess impression material.



2. Mount two supporting pins and pour counter impression.



3. Mount two pins per prep and supporting pins on adjacent teeth.

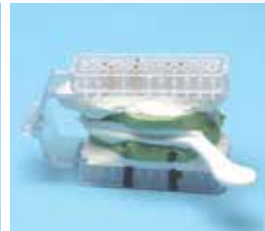

4. Pour master impression.



5. Remove impression tray and section dies.

Ordering Info

0069	W.O.W. Intro Kit
0065	W.O.W. Full Arch
0065AD	W.O.W.2 Full Arch
0072	W.O.W. Posterior
0072AD	W.O.W.2 Posterior
0067	Dowel pins (300)
0020	Die Maker

Intro kit contains 20 full arch, 30 posterior and 300 dowel pins.
All other W.O.W. kits contain 50 articulators and 300 dowel pins.
Die maker comes complete with 50 posterior articulators.

Provide a low cost, professionally presented case to your customers with the W.O.W. Articulators!

CE

EC REP
MDSS GmbH
Schiffgraben 41
30175 Hannover, Germany

#1

TSZTECH™

The P.B.S.® Brazed Advantage

Faster, Cooler Cutting and Longer Lasting!

A New Diamond Technology for Today's Advanced Materials - Lithium Disilicate and Zirconia

Highest rated all-ceramic cutting for today's new generation of tougher ceramic materials.

Those materials include but are not limited to IPS e.max®, BruxZir®, Lava™ Plus, Cercon®, ZirLux® FC, Obsidian™, IPS Empress, Procera® and Lava™ Ultimate.

TSZtech Contraindications – Will not cut or separate metal.



Laboratory Diamonds
Available one per pack:

2015895 **TSZtech™ Diamond Assorted 5-pack**
(260.8Z, 285.5Z, 703.8KRZ, 125Z, 770.8Z)

5-pack Reorder Item Numbers:



260.8Z
2015893
Minor/Major
Diameter
0.3 / 1.3



285.5Z
2015891
Minor/Major
Diameter
2.2



703.8KRZ
2015894
Minor/Major
Diameter
1.1 / 1.8



125Z
2015892
Minor/Major
Diameter
1.7



770.8Z
2015890
Minor/Major
Diameter
1.2 / 1.7



"L" - 1.5mm Longer Shank
Available in 5-packs:



L260.8Z
2015883
Minor/Major
Diameter
0.3 / 1.3



L125Z
2015882
Minor/Major
Diameter
1.7



L770.8Z
2015880
Minor/Major
Diameter
1.2 / 1.7



X590.12Z
2012847
Minor/Major
Diameter
1.2



X769.10Z
2012845
Minor/Major
Diameter
0.7 / 1.8

Two Striper®
An abrasive technology product

Distributed by
premier

Premier® Dental Products Company • 1710 Romano Drive, Plymouth Meeting, PA 19462, USA

Phone 888-670-6100 • 610-239-6053 • Fax 610-239-6171 • www.premusa.com • E-mail: dentalinfo@premusa.com

© Premier is a registered trademark of Premier Dental Products Company. / ®Two Striper is a registered trademark and TSZtech is a trademark of the product manufacturer, Abrasive Technology, Inc. / U.S. Patent 5, 511, 718.
© IPS e.max, BruxZir, Cercon, ZirLux FC, IPS Empress, Procera are not registered trademarks and™ Lava Plus, Lava Ultimate and Obsidian are not trademarks of Premier Dental Products Company.

8621290 051359 Rev2 NKP 2M