

**Reliably the Best for  
Accurate Impressions**



**19x  
winner!**

## Triple Tray®

### Dual-Arch Impression Trays

Premier Triple Tray® is the most trusted name for accurate dual-arch impressions – proven over 185 Million impressions and counting. The decision is simple. The superior performance of Triple Tray makes it the go-to tray for accurate impressions the first time, every time.

**3→1**

Three techniques in one  
*Saves time, money and materials*



Adaptable, ultra-thin webbing  
*Provides accurate bite registration and low water sorption<sup>2</sup>, reducing the possibility of distortion*



Disposable  
*No cleaning involved, aseptic*



Seven configurations  
*Designed to address all clinical cases*



1006283 Wide Body Posterior (48 trays/box)

1006243 Sideless Posterior (48 trays/box)

1006203 Posterior (48 trays/box)

1006253 Extended Quadrant (35 trays/box)

1006213 Anterior (35 trays/box)

1006273 Three-Quarter Arch (35 trays/box)

1006223 Quadrant (40 trays/box)

1006233 Assorted (5 each of the styles above)

**premier®** Inspired solutions  
for daily dentistry®

Premier® Dental Products Company • 888-670-6100 • premierdentalco.com

Ensure accuracy and save time

## T-LOC™ Triple Tray®

Adhesive-Free Dual-Arch Impression Trays

- ✓ **Patented adhesive-free design** - Retention tracks mechanically lock in impression material
- ✓ **Easier handling** - Large, ergonomic handles simplify loading and transferring the tray
- ✓ **Accurate bite** - Thin, loose webbing maximizes intercuspation of the teeth to ensure an accurate bite
- ✓ **Ensure accuracy and save time** - Low sides, thin posterior bar and wide impression areas avoid impingement

1006274 Three-Quarter Arch  
(22 trays/box)

1006204 Posterior  
(35 trays/box)

1006214 Anterior  
(28 trays/box)

U.S. Patent #6,749,428



### Patented Design

Retention tracks mechanically lock in material - from all sides - eliminating need for adhesive.



## Alfa™ Triple Tray®

Full-Arch Alginate Impressions

- ✓ **Revolutionizing full-arch alginate impression making**

Alfa™ Triple Tray® improves patient comfort while helping to prevent gagging. Alfa cuts chair time in half and is easy to use for most full-arch alginate impressions. Use for custom trays, whitening trays, night guards, athletic guards, study models, and as a matrix for temporaries.

1006313 Small (purple)  
(24 trays/box)

1006323 Medium (white)  
(24 trays/box)

1006333 Large (green)  
(24 trays/box)

1006343 Assorted  
(5 small, 14 medium, 5 large)

U.S. Patent #6,749,428




**Retention Tracks**  
eliminate adhesives

**Purchase Point**  
on each side for  
easy removal

**Thin, Loose  
Webbing**  
for complete  
intercuspation

**Palate-Free Design**  
prevents gagging

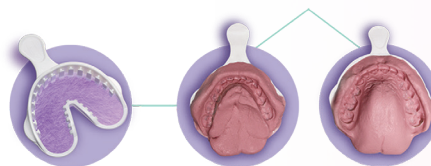

One Tray - Two Impressions

**premier®** Inspired solutions  
for daily dentistry®



Premier® Dental Products Company • 1710 Romano Drive, Plymouth Meeting, PA 19462, USA  
Phone 888-670-6100 • 610-239-6000 • Fax 610-239-6171 • premierdentalco.com • E-mail: info@premierdentalco.com

\* Premier and Triple Tray are registered trademarks and Alfa and T-LOC are trademarks of Premier Dental Products Company. / ©2023 Premier Dental Products Company.  
1. 2003, 2004, 2006, 2007, 2008, 2009, 2010, 2011, 2012, 2013, 2014, 2015, 2016, 2017, 2018, 2019, 2020, 2021, 2022 Townie Choice Award, DentalTown Magazine. / 2. Independent Research Organization, July 2000.

8698741 0623034 Rev7 PDF