

New!

Introducing! Traxodent® Gels



Traxodent® Gel AlCl

25% Aluminum Chloride Hemostatic System

***Clear Control.
Proven Performance.***

Traxodent® Gel AlCl delivers fast, reliable hemostasis and sulcular fluid control. The clear non-staining formula rinses clean and does not leave a residue that can interfere with impressions, bonding and cementation procedures. Use Traxodent Gel AlCl for crown & bridge, composite restorations, impressions, or anytime hemostasis is required.



Traxodent® Gel FS

20% Ferric Sulfate Hemostatic System

***Targeted Hemostasis with
the Strength of Ferric Sulfate.***

Traxodent® Gel FS is a powerful 20% ferric sulfate hemostatic gel formulated for fast, reliable bleeding control during restorative and prosthetic procedures. The gel's optimal handling makes it easy to place exactly where it's needed and helps prevent gel spreading to nearby tissue.



premier® Inspired solutions
for daily dentistry®

Premier® Dental Products Company • 888-670-6100 • premierdentalco.com

©2025 Premier Dental | All Rights Reserved

Traxodent® Gel AlCl

25% Aluminum
Chloride
Hemostatic System



NEW!

Traxodent® Gel FS

20% Ferric Sulfate
Hemostatic System



***Clinically proven results you
can count on.***

- ✓ 25% aluminum chloride for effective hemostasis & fluid control
- ✓ Innovative formula boosts gel penetration for faster, more effective results
- ✓ Excellent visibility during placement & rinse
- ✓ Syringe delivery for easy, accurate application

9007101 Traxodent® Gel AlCl Starter Kit

Contents: 30cc (30mL) Filled Bulk Syringe,
20 Brush Tips, 20 x 1.2ml Syringes,
20 x 1.2ml Syringe Caps, Docking Port

9007105 Traxodent® Gel Brush Tips

Contents: 100 Brush Tips

***Controlled performance
where it counts.***

- ✓ 20% ferric sulfate gel with optimized handling
- ✓ Rapid hemostasis & sulcular fluid control for procedural efficiency
- ✓ Excellent visibility during placement & rinse
- ✓ Syringe delivery for easy, accurate application


9007106 Traxodent® Gel FS Starter Kit

Contents: 30cc (30mL) Filled Bulk Syringe,
20 Brush Tips, 20 x 1.2ml Syringes,
20 x 1.2ml Syringe Caps, Docking Port


9007105 Traxodent® Gel Brush Tips

Contents: 100 Brush Tips


Easy-to-Use Technique

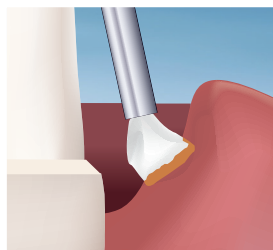

1 Work Traxodent Gel firmly into the sulcus with the brush tip until bleeding is controlled.

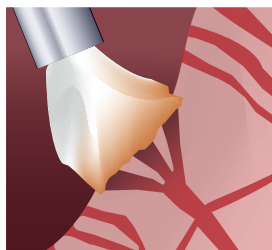

2 The brush tip enables gentle pressure, helping Traxodent Gel penetrate cut capillaries.

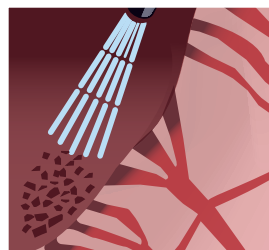

3 Rinse with air/water to remove the coagulum. If bleeding continues, reapply.

Available through your authorized dealer.

premier® Inspired solutions
for daily dentistry®



Premier® Dental Products Company • 1710 Romano Drive, Plymouth Meeting, PA 19462, USA
Phone 888-670-6100 • 610-239-6000 • Fax 610-239-6171 • premierdentalco.com • E-mail: info@premierdentalco.com

Lit# 8610438 0625067 Rev0 NKPSM