


New!

Precise Cure

Technology Inspires
Clinical Confidence

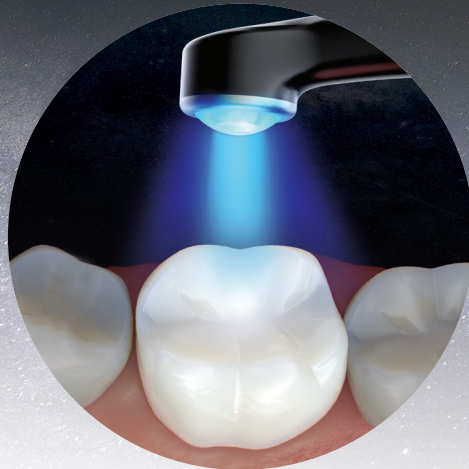

Lumira™ LED Curing Light

Designed for today's clinicians, the Lumira™ LED curing light delivers consistent performance and durability. Lumira's lightweight, all-metal construction withstands the rigors of daily dental use while reducing hand fatigue. Its innovative low-profile head rotates 360° for effortless access and optimal illumination throughout the oral cavity. Lumira combines advanced optical technology with clinical versatility, offering a broad wavelength range to effectively polymerize CQ products, a maximum output of up to 3000 mW/cm² for rapid curing, and a high-intensity collimated beam for uniform light distribution and increased depth of cure.

With seven curing modes – including Normal, High, Turbo, Ortho, Soft Pulse and Check (Transillumination) – the Lumira LED curing light provides a comprehensive solution for delivering predictable, reliable cures, even in challenging deep restorations.

premier® Inspired solutions
for daily dentistry®

Premier® Dental Products Company • 888-670-6100 • premierdentalco.com

©2026 Premier Dental | All Rights Reserved

Lumira™

LED Curing Light



Wave Range

- 440–490 nm
- Polymerizes CQ products



Intensity

- Max irradiance 3000 (mW/cm²)
- Matches premium brands



Collimated Beam

- Less light loss¹
- Reliable depth of cure



7 Clinical Modes

- Align with restorative workflows
- Normal, High, Turbo, Pulse, Soft, Ortho & Check



Streamline Design

- Durable all metal design
- Ergonomic & lightweight (108g)



Ergonomic & Easy to Use

Simple 2 button design with OLED screen displays mode, battery life, light intensity, and time clearly

Charging Base with Built-in Light Meter

Ensure accurate light intensity for efficient polymerization




3006100 Lumira™ LED Curing Light System

Contents:

- 1 x Lumira™ Curing Light Unit
- 1 x LED Light Head
- 1 x Charging Station
- 1 x Adapter
- 1 x Light Shield
- 1 x Point Cure Lens

Available through your authorized dealer.


Less Light Loss¹

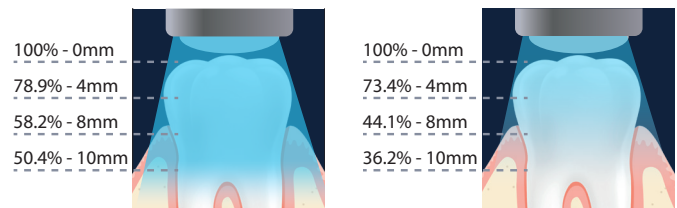


Lumira™

Leading Brand

Low Profile, 360° Rotatable Head for Easier Access

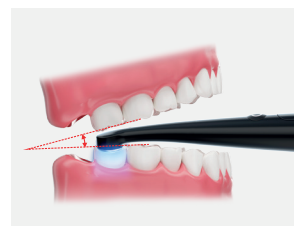

Lumira™

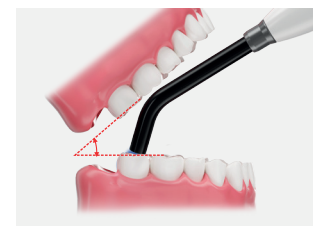

Leading Brand

premier® Inspired solutions for daily dentistry®



Premier® Dental Products Company • 1710 Romano Drive, Plymouth Meeting, PA 19462, USA
 Phone 888-670-6100 • 610-239-6000 • Fax 610-239-6171 • premierdentalco.com • E-mail: info@premierdentalco.com

1. Data on file.